

## USEFUL INFORMATION



- Hair follicles that have been deactivated by the Tria Laser 4X will not grow back.
- Laser hair removal does not cause hair to fall out immediately. After your first treatment, there will be no visible change to the hair.
- Some users may feel the sensation of a rubber band snapping on the skin. Others will feel a warm, prickling sensation. As the amount of hair decreases, the sensation will also decrease.
- In a clinical study, higher levels were found to be 2-3 times more effective. So for best results, we recommend you use the highest setting that you can tolerate.
- Treat every 2 weeks for up to 3 months or until you reach the desired results.

Clinical studies on file

## CHECK YOUR SKIN TONE

The Tria Laser 4X is intended only for light to medium skin tones and medium to dark hair. It does not work on white, grey, blonde or red hair. If your skin tone matches any of the skin tones labeled "Not Safe" for treatment, then DO NOT USE the Tria Laser 4X.

### SKIN TONE CHART

Intended Skin Tones:



White    Ivory    Beige    Light Brown

Not Safe:



Medium Brown    Dark Brown or Black

**FOR USE ON** Legs • Arms • Bikini Line • Underarms • Feet • Hands • Chest • Back • Stomach • Women's Facial Hair (cheek line downward)

### TRIA BEAUTY CUSTOMER CARE

United States  
877-321-8742

customercare@triabeauty.com  
triabeauty.com

### United Kingdom/Europe

0800 0484 322


info@triabeauty.co.uk  
triabeauty.co.uk

tria.

**HAIR REMOVAL LASER 4X**

Quick Start Guide

## GET TO KNOW YOUR TRIA LASER 4X

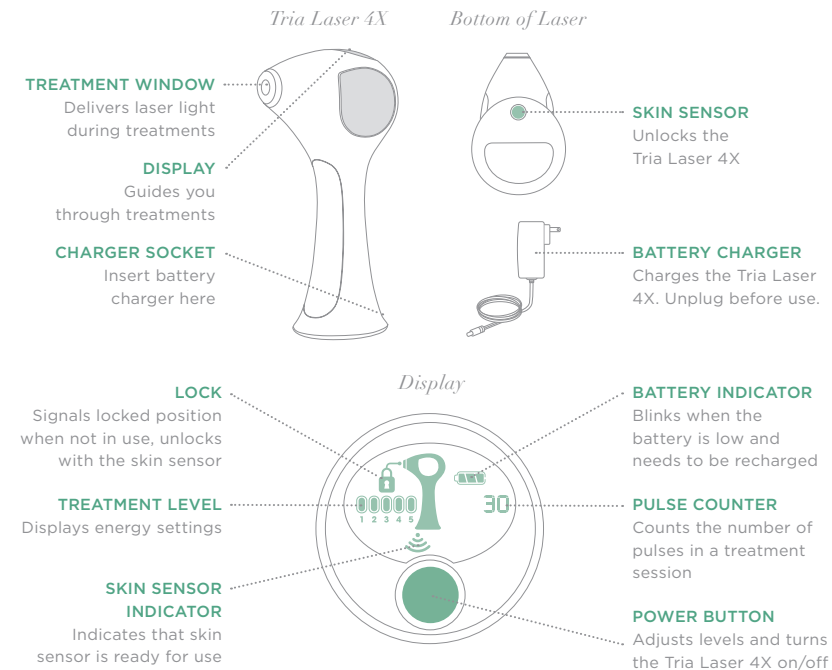

## BEFORE FIRST TREATMENT

For detailed directions and important safety warnings, please read the complete *Instructions for Use* guide before use.

### HOW TO TREAT

Before every treatment, shave, clean and dry the area you wish to treat. Before your first full treatment on any new area, perform a patch test in that area. Repeat treatments every 2 weeks for up to 3 months or until you're satisfied with results.


#### 1 TURN ON DEVICE

Unplug the Tria Laser 4X from the battery charger and press the power button.



#### 2 UNLOCK

Unlock the device by holding the Skin Sensor against the skin in the area you wish to treat. If the Tria Laser 4X is safe for your skin tone, a tone will sound and the Tria Laser 4X will unlock (if device doesn't unlock, see *Instructions for Use* for more information).

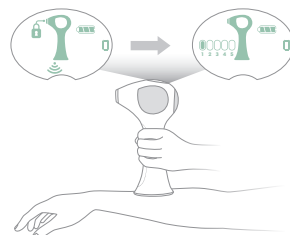


**CHARGE:** Charge the Tria Laser 4X until the battery light is full (~2 hours). Unplug before use.

**ACTIVATE:** A one-time, 5-minute activation is required before your first treatment. In the US: visit [tria beauty.com/activate](http://tria beauty.com/activate) or call Customer Care at 1-877-321-8742. In the UK: visit [tria beauty.co.uk/activate](http://tria beauty.co.uk/activate) or call Customer Services at 0800 0484 322.


#### 3 SET THE TREATMENT LEVEL

Press the power button to adjust levels. For best results, use the highest setting that you can tolerate.

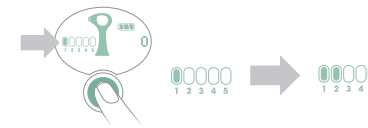


#### 4 TREAT. HOLD FOR 2 BEEPS.

- Touch and hold the treatment window to the skin until you hear two beeps.
- After the two beeps have sounded, lift the laser and place it back on the skin so that it overlaps the previous spot (see overlapping pattern diagram).

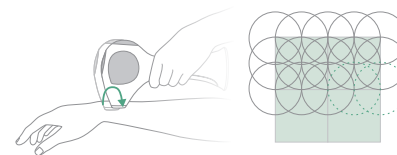


#### 5 KEEP GOING TO COVER THE AREA

Use the Tria Laser 4X Pulse Counter as a guide to how many pulses you may need to cover each area per treatment.

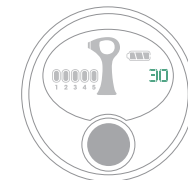

## PULSE COUNTER GUIDE

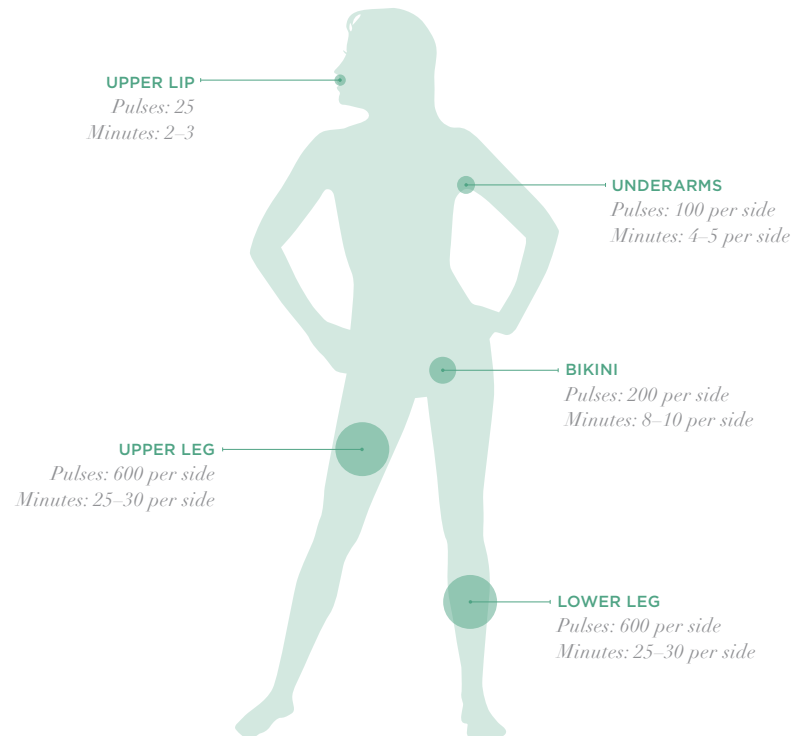

Pulse counts and times are only approximations—your actual treatment area and number of pulses may vary. Always follow the treatment instructions and warnings provided in the *Instructions for Use* booklet. When treating larger areas like legs and back, you may need more than 600 pulses to cover the area.